

Cable No: FP7

Length: 34 cm

Cable Name: cable (FP7) of Assembly of Red wire-one side Golden lug & other side Batt contact clip

Cable Type: 1 Sq mm

Note:


After crimping battery contact clip, put solder on it. After this use heat shrink cables sleeves & then push the plastic cap on clip.

Placement of the cable

Golden lug – U plate of Sternum handle (Backside of the front panel, on the Sternum handle)

Battery contact clip – On PCB 12026 ident (CON 5)

Diagram:

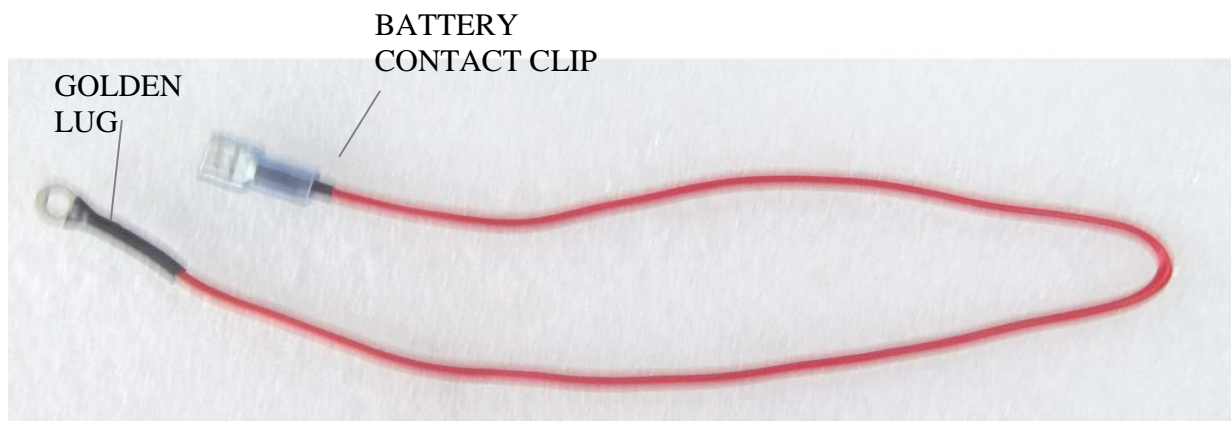

TOTAL LENGTH: FP7-34CM