

Cable SP4

Note: Use Heat shrink sleeve on both sides to all cables.
 Cable Length is without connector.
 Color code given is recommended color code.
 Twist the cables.


**FP - Front Plate, BP - Back Plate, PS - Power Supply, DC – Data Cable
 SP – Side Plate**

Cable : -SP4 NIBP Cable Length: 30 cm

(Lifeguard NIBP Module to 6 pin Relimate on PCB2K12032)

Specification of Connector 1: 6 Pin JST [F] 2.54 mm on Lifeguard NIBP PCB

Specification of Connector 2: 6 Pin Relimate [F] 2.54mm on J5 PCB 2K12032

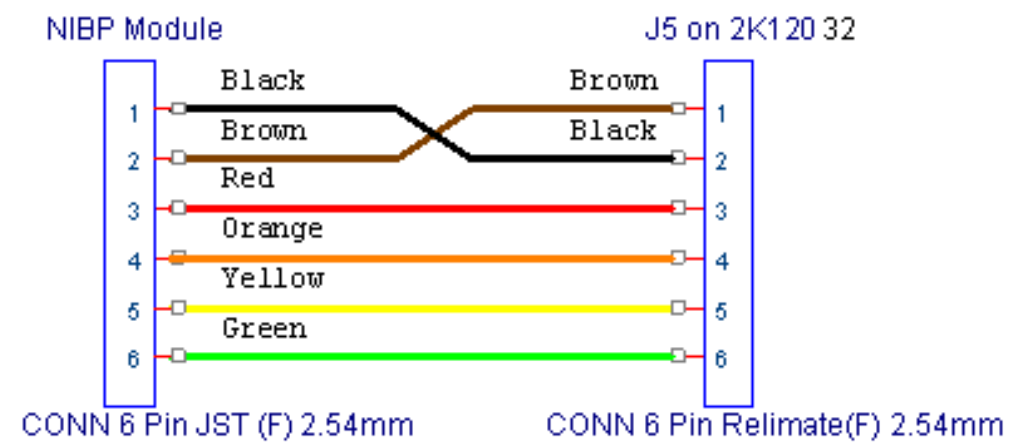

Connector 1	Wire Color (4 core shielded cable)	Connector 2
1	Black	2
2	Brown	1
3	Red	3
4	Orange	4
5	Yellow	5
6	Green	6