



DEDICA

Dedicated Stress Test System

Key Features

- ✓ 17" colour TFT display with Touch screen operations
- ✓ Offers Rock steady baseline at heart rate more than 210
- ✓ Offers full frequency response from 0.05 to 150 Hz. as stipulated by American Heart Association.
- ✓ Offers accurate computation of ST levels, assuring you of excellent correlation with angiography.
- ✓ Offers clean ECG without muscle tremor or 50 cycle noise.
- ✓ ST integral of all 12 leads

Technical Specifications

General

Length	: 2055 Mm.
Width	: 750 Mm.
Walking Area	: 1440 X 500 Mm.
Weight	: 100 Kg.
Belt Speed	: Smooth, 0 - 15 Kmph.
Incline Range	: 0 - 22 %
Patient Weight	: Upto 250 Kg.
Power	: 230v, 50hz, 15a
In-house Manufacturing	: Yes



Features

- ✓ During acquisition when in auto mode, the enlarged median having maximum ST depression is automatically selected and displayed.
- ✓ Enlarged median can be selected manually as well.
- ✓ During acquisition, medians can be refreshed at one second interval from 1 to 8 seconds.
- ✓ Data can be reviewed beat to beat with full disclosure for all 12 leads of the entire test.
- ✓ There is virtually no constraint on the number of cases which can be stored.

Features of TreadStar XP

- ✓ Our treadmills accelerate slowly from a standstill and gradually gather speed. They have high torque at low speed.
- ✓ Wide and long belt for stable and comfortable stride. One-inch wooden board for extended life.
- ✓ Anti-skid side steps.
- ✓ Powerful DC motor with 6 HP peak power capabilities can ensure trouble-free running under most heavy loads.
- ✓ AC motor driven integrated front incline system with built in linear actuator (1100 Lbs thrust) incorporates feedback for correct incline positioning.
- ✓ Smooth, noise free elevation.



[Click Here For Technical Specification](#)

*Due to our continuous product development program, the specifications may change without prior notice.


12 Lead mode with Superimposed Medians



4 Lead mode with ST Integral



3 Lead mode

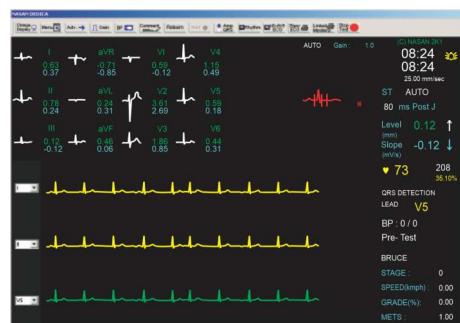

NASAN MEDICAL ELECTRONICS PVT. LTD.

C Wing - Shubham Heights, Off Mumbai - Bangalore Highway Service Road, Warje, Pune - 411058

Phone : +91 8237032389 | E-mail : sales@nasanmedical.com

www.nasanmedical.com