

Cable No: FP3


Cable name:-Assembly of Encoder for Defib

Specification of Connector 1:- 2.54mm SHS (F) 6 pin Relimate connector J7 on PCB 2k13002

Specification of Connector 2:- Mechanical Encoder

Length -30 cm

Cable Diagram:

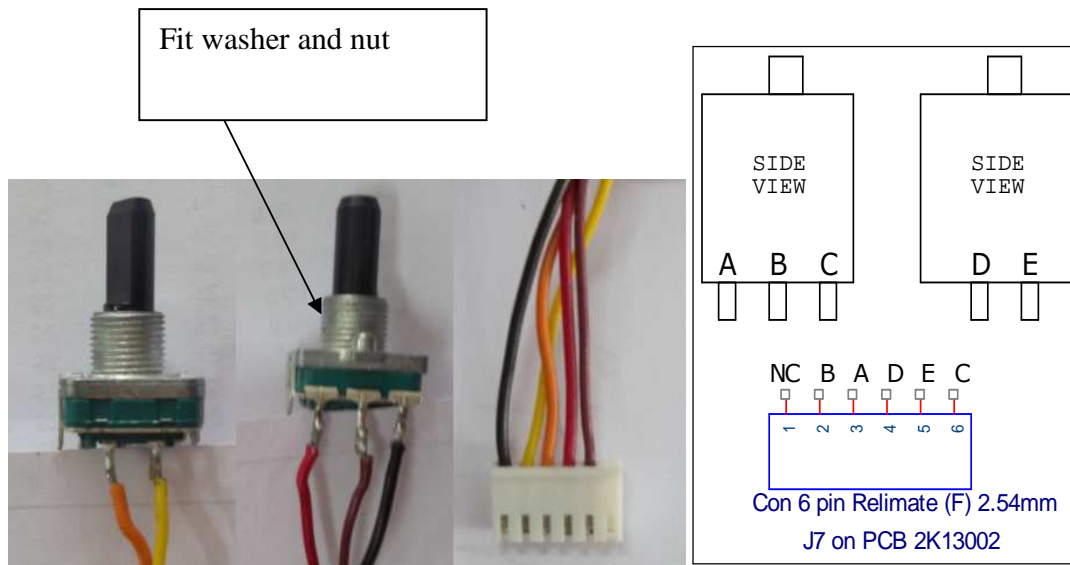

Connection details:

6 Pin SHS Relimate [F] Connector 2.54 mm J7 on 2K13002	7/36 PVC Wire Color Code	Mechanical Encoder
1	Black	NC
2	Brown	B
3	Red	A
4	orange	D
5	yellow	E
6	Green	C

After connection fit plane washer and nut which comes along with encoder.