

Cable SP 2.2


Note: Use Heat shrink sleeve on both sides to all cables.
 Cable Length is without connector.
 Color code given is recommended color code.
 Twist the cables.

FP - Front Plate, BP - Back Plate, PS - Power Supply, DC – Data Cable
SP – Side Plate

Cable : SP 2.2- SPO2 Connection Length: 30 cm

Specification of Connector 1: 4 Pin JST [F] 2.54mm on SpO2 PCB 14015

Specification of Connector 2: 4 Pin Relimate [F] on CN24 PCB 2K9013

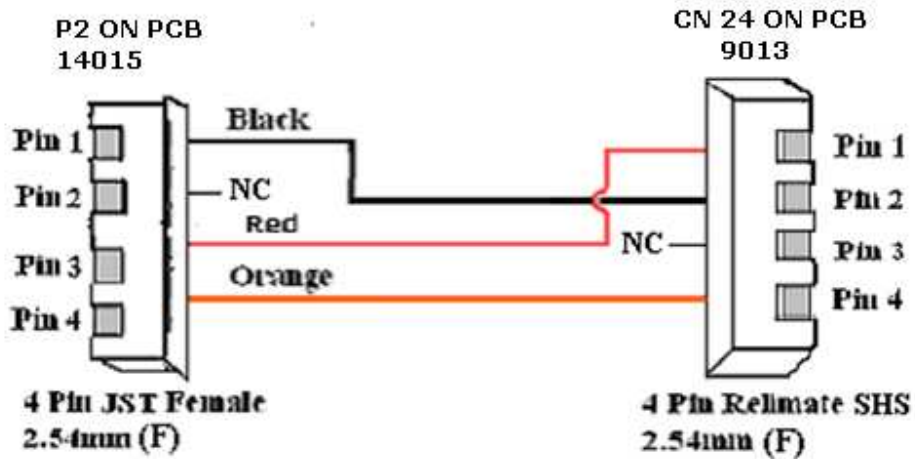

Connector 1	Wire Color (4 core shielded cable)	Connector 2
1	Black	2
2		NC
3	Red	1
4	Orange	4